



COUCHES AND STRETCHERS

GYNECOLOGY AND UROLOGY COUCH ADJUSTABLE HEIGHT

Table specially designed for gynecology and urology examinations.

Highest quality upholstery, very easy to clean. Fire retardant and antibacterial, resistant to urine, blood and disinfectant products, specially recommended for its use in hospitals.

The armchair is equipped with 3 motors, very easy to operate by the patient and the medical staff.

Backrest and footrest can be moved both independently or simultaneously till horizontal position (1950 mm)

Adjustable height range from 550 to 900 mm.

Trendelenburg position 12°.

Backrest includes an adjustable pillow for patient's comfort.

Seat specially designed for gynecology with a stainless steel waste tray under seat.

Includes anatomical upholstered legrests, adjustable in height and position.

Four ø100 mm swivel castors. All of them have break.

Structure protected by an ABS plastic shell.

Paper roll holder is included.

Structure manufactured in enamelled epoxy steel.

Safe working load: 175 Kg.

· Bed dimensions: 1950 x 600 mm

· Backrest: 900 x 600 mm

· Seat: 500 x 600 mm

· Footrest: 500 x 600 mm



H-301 .. 1950 - 900 x 670 x 1400 mm .. 75 kg

GYNECOLOGICAL COUCH FIXED HEIGHT

Articulated backrest and footrest from by mechanical ratchet from 90° to 180°.

Trendelenburg position up to 45°, with security mechanism.

Upholstered in vinyl, colour to choose, fireproof, with high density sponge foam inside.

Set of anatomical leg rests, adjustable in height and direction.

Moulded seat with waste tray underneath.

Legs with non-slipping ends.

Structure made in steel tube with dimensions $\varnothing 30 \times 1,5$ mm thickness, with reinforcing cross-piece
Available finishings: stainless steel 18/10 or epoxy coating.
Paper roll holder included as standard.
Safe working load: 150 Kg

· Bed dimensions: 1550 x 500 mm

· Backrest: 680 x 500 mm

· Seat: 380 x 500 mm

· Footrest: 490 x 500 mm

Trendelenburg position 45°

H·02



EXAMINATION COUCH FIXED HEIGHT

Backrest and footrest articulated in three positions by mechanical ratchet.

Upholstered in vinyl, colour to choose, fireproof, with high density sponge foam inside.

Legs with non-slipping ends.

Structure made in steel tube with dimensions $\varnothing 30 \times 1,5$ mm thickness, with reinforcing cross-piece

Available finishings: stainless steel 18/10, chrome or epoxy coating.

Available accessories: paper roll holder and paper cutter.

· Bed dimensions: 1800 x 550 mm

· Backrest: 500 x 550 mm

· Seat: 750 x 550 mm

· Footrest: 550 x 550 mm

H·08

three-sections model



H·02 .. 1550 - 1050 x 650 x 850 mm .. 32 kg

H·08 .. 1800 - 1250 x 550 x 750 mm .. 23 kg



COUCHES AND STRETCHERS

EXAMINATION COUCHES

FIXED HEIGHT

Articulated headrest in three positions by mechanical ratchet (45°).

Upholstered in vinyl, colour to choose, fireproof, with high density sponge foam inside.

Legs with non-slipping ends.

Structure made in steel tube with dimensions $\varnothing 30 \times 1,5$ mm thickness, with reinforcing cross-piece

Available finishings: stainless steel 18/10, chrome or epoxy coating.

See available accessories on page 9.

two-sections model

H-03



optional accessories
" see page 9 "



H-03 EXAMINATION TABLES FOR ADULTS

1800 x 600 x 700 mm ∴ 24 kg

1900 x 600 x 700 mm ∴ 25 kg

1900 x 600 x 800 mm ∴ 27 kg

1900 x 700 x 700 mm ∴ 30 kg

1900 x 700 x 800 mm ∴ 33 kg

H-07 EXAMINATION TABLE FOR CHILDREN

1500 x 500 x 800 mm ∴ 18 kg

BABY CHANGING STATION

Ergonomic platform with smooth and rounded shape without sharp edges.

Wide and deep surface with adjustable safety belt for safe use and easy handling.

Convenient opening and closing with one hand.

Folding system with safety lock when not in use.

Possibility of hanging toilet bag on each side.

Model to be installed on the wall with four fixing bolts for maximum safety.

Manufactured in reinforced polymers with high impact resistance, with phytosanitary treatment.

Safe working load: 15 kg

DIMENSIONS

Open: 820 x 540 x 480 mm (length x width x height)

Closed: 820 x 115 x 480 mm (length x width x height)



H·09



wall mounted

PAEDIATRIC EXAMINATION TABLE

Upholstered in vinyl, colour to choose, fireproof, with high density sponge foam inside.

Legs with non-slipping ends.

Lower shelf made of laminated wood.

Available finishings: stainless steel 18/10 or epoxy coating.

bed with lateral protection

H·10



H·10 .. 1200 x 600 x 850 mm .. 30 kg



COUCHES AND STRETCHERS

EXAMINATION COUCHES

ADJUSTABLE HEIGHT

“Scissors” design for vertical lifting, without lateral movement to minimise space and provide maximum stability throughout all profiling positions.

Hydraulic or electric system for raising and lowering.

Height range: 480 - 1050 mm

Articulated backrest up to +80° bay gas spring.

Backrest with breathing hole and cap.

Upholstered in vinyl, colour to choose, fireproof, with high density sponge foam inside.

Retractable wheeling system.

Structure made in epoxy coated steel.

Paper roll holder is included as standard.

- Bed dimensions: 1900 x 700 mm
- Backrest: 700 x 700 mm
- Seat: 1200 x 700 mm
- Safe working load: 200 kg

H·04



retractable castors


H·04 ∴ 1900 x 700 x 480 - 1050 mm ∴ 65 kg

ACCESSORIES FOR EXAMINATION COUCHES



H-901

paper roll holder in ABS

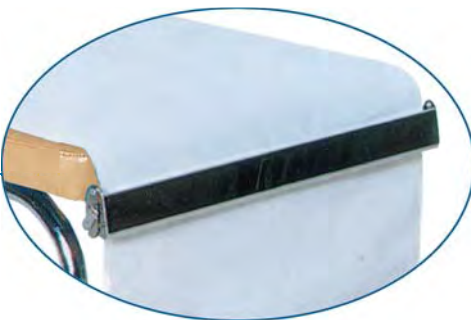


H-902

paper cutter in stainless steel




H-905

set of 4 castors

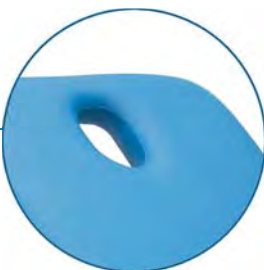

H-903


breathing hole

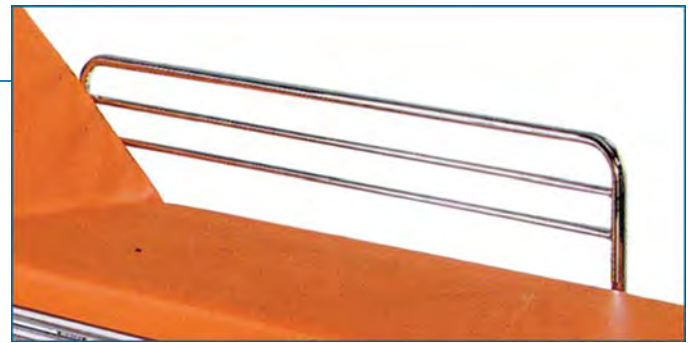

H-931

security rails

supporting table

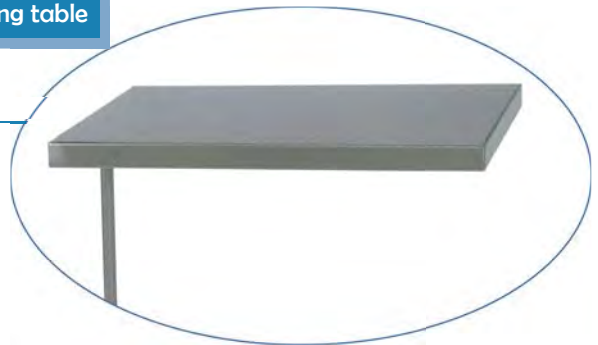

H-938

armrest



H-939



H-935

infusion stand



H-951

set of anatomical leg holders

H-901 .. Paper roll holder

H-902 .. Paper cutter

H-903 .. Breathing hole with cap

H-905 .. Set of 4 castors

H-931 .. Stainless steel security rails

H-935 .. Stainless steel infusion stand

H-938 .. Supporting table

H-939 .. Articulated armrest

H-951 .. Set of anatomical leg holders



COUCHES AND STRETCHERS

TREATMENT COUCHES

ADJUSTABLE HEIGHT

Hydraulic or electric system for raising and lowering.

Height range: 550 - 900 mm

Articulated backrest +80° and -30° by gas spring.

Headrest with breathing hole and cap.

Upholstered in vinyl, colour to choose, fireproof, with high density sponge foam inside.

Retractable wheeling system.

Structure made in epoxy coated steel.

Paper roll holder is included as standard.

.. Bed dimensions: 1900 x 700 mm

.. Headrest: 700 x 700 mm

.. Body section: 1200 x 700 mm

.. Safe working load: 200 Kg

H-05

hydraulic

double pedal on each side

two sections

electrical

retractable castors

H-05 .. 1900 x 700 x 550 - 900 mm .. 65 kg

TREATMENT COUCHES

ADJUSTABLE HEIGHT

Hydraulic or electric system for raising and lowering.
 Height range: 550 - 900 mm.
 Articulated backrest +70° and -15° by gas spring.
 Articulated footrest +30° and -30° by gas spring.
 Headrest with breathing hole and cap.
 Upholstered in vinyl, colour to choose, fireproof, with high density sponge foam inside.
 Retractable wheeling system.
 Structure made in epoxy coated steel.
 Paper roll holder is included as standard.

- Bed dimensions: 1900 x 700 mm
- Backrest: 700 x 700 mm
- Seat: 600 x 700 mm
- Footrest: 600 x 700 mm
- Safe working load: 200 Kg



hydraulic

H-06

double pedal on each side

three sections

electrical

retractable castors

H-06 .. 1900 x 700 x 550 - 900 mm .. 70 kg



COUCHES AND STRETCHERS

TREATMENT COUCHES ADJUSTABLE HEIGHT

Hydraulic system for raising and lowering actuated by pedals (one at each side).

Height range: 550 - 900 mm

Trendelenburg (20°) and reverse trendelenburg (30°) actuated by progressive gas spring with lateral levers (one at each side)

Articulated backrest by gas spring (80°), with breathing hole and cap.

Upholstered in vinyl, colour to choose, fireproof, with high density sponge foam inside.

Retractable wheeling system.

Side rails for placing accessories (see page 9)

Lower shelf for patient's belongings and structure shell.

Structure made in epoxy coated steel.

Paper roll holder is included as standard.

· Bed dimensions: 2000 x 700 mm

· Backrest: 700 x 700 mm

· Asiento: 1300 x 700 mm

· Safe working load: 200 kg

Trendelenburg position



Reverse Trendelenburg position



H·305

lateral rails for accessories



retractable castors

hydraulic

H·305 .. 2000 x 700 x 550 - 900 mm .. 90 kg

TRANSPORT STRETCHER ADJUSTABLE HEIGHT

Two sections bed base with articulated headrest (60°) actuated by gas spring.

Upholstered in vinyl, colour to choose, fireproof, with high density sponge foam inside.

Hydraulic system for raising and lowering, actuated by pedal.

Reinforced side rails, with quick and accurate folding system.

I.V. pole and oxygen bottle holder included.

Handle for easy movement.

ø150 mm castors with directional system and centralized brake. Structure made in epoxy coated steel.

- Bed dimensions: 1900 x 600 mm
- Height range adjustment: 600 - 1200 mm
- Safe working load: 200 kg

gas spring assisted backrest

reinforced folding side rails

H-33

hydraulic

ø150 mm castors

directional system and centralized brake



H-33 .. 2000 x 700 x 600 - 1200 mm .. 72 kg



COUCHES AND STRETCHERS

TRANSPORT STRETCHERS

FIXED HEIGHT

Handle for easy movement.
 Upholstered in vinyl, colour to choose, fireproof, with high density sponge foam inside.
 Reinforced side rails, with quick and accurate folding system.
 I.V. pole and oxygen bottle holder included.
 Structure made in stainless steel 18/10.

H-35

Articulated headrest (45°) by mechanical ratchet in three positions.
 ø125 mm castors, two of them with brake, with rubber bumper rings.
 .. Bed dimensions: 1900 x 550 mm

H-34

Articulated headrest (60°) by gas spring with actioning lever.
 ø150 mm castors, two of them with brake, with rubber bumper rings.
 .. Bed dimensions: 1900 x 550 mm

H-35

manual backrest tilting



ø125 mm castors

H-34

reinforced folding side rails

gas spring assisted backrest



ø150 mm castors

H-34 .. 2000 x 700 x 850 mm .. 37 kg

H-35 .. 2000 x 700 x 850 mm .. 35 kg

EMERGENCY STRETCHER

Excellent capacity of transport and easy operation.
 X-ray permeable bed with detachable mattress with dimensions 1950 x 650 mm, with 100 mm thick sponge foam core upholstered in fireproof vinyl with antibacterial, anti-mite and antifungal treatment, resistant to abrasion and stains.
 Articulated backrest (70°) assisted by gas spring, with lever at each side.
 Height adjustable by hydraulic pump actuated by double pedal.
 Trendelenburg (20°) and reverse trendelenburg (30°) system, assisted by progressive gas spring, with levers at each side.
 Folding rails with easy fast and safe handling with security block/unblock system.
 Lateral handles for easy movement.
 Included accessories: object bearing tray, two I.V. poles and oxygen bottle holder.
 Other optional accessories: Monitor shelf, instrument tray, and spare I.V. pole.

Set of four ø200 mm swivel castors with directional system and centralized brake.
 Directional supporting fifth wheel.
 Perimetrical protection bumpers.
 Structure made in epoxy coated steel. Accessories made in chromed and stainless steel.

- Bed dimensions: 1950 x 650 mm
- Height range: 650 - 1000 mm
- Safe working load: 200 kg



Trendelenburg and Reverse Trend. position

transversely folding rails



H-38

fifth wheel

ø200 mm castors with directional system and centralized brake

H-38 .. 2200 x 860 x 650 - 1000 mm .. 115 kg